



**MILLENNIUM
SURGICAL**


SURGICALINSTRUMENTS.COM

UNMATCHED SERVICE · UNLIMITED INSTRUMENTS

822 Montgomery Ave #205
Narberth, PA 19072
(800) 600-0428


Fallopian Ring Applicator

Female sterilization by means of Silastic Ring Applicator is a clinically well-established procedure based on decades of experience. Our instruments are completely reusable and fully autoclavable. They are manufactured of the highest grade materials by experienced and skilled craftsmen.

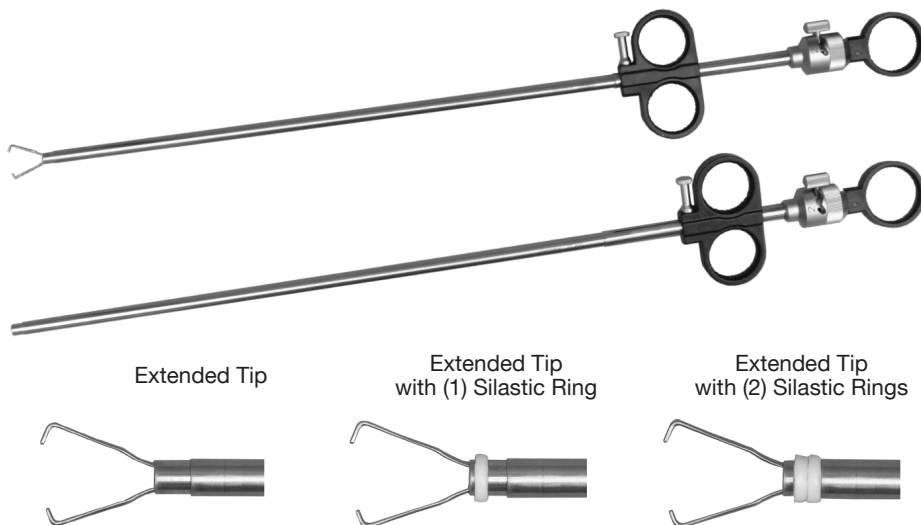

The Silastic Ring Applicators are made for single or dual ring application with one single insertion. The instruments are available in three different lengths for use through the operating scope (single puncture approach) or through trocar sleeve in conjunction with a laparoscope (multi-puncture approach).

LIST NO. 52105-500	6MM Ø, 40CM SILASTIC RING APPLICATOR- FOR SINGLE AND DUAL RING APPLICATION; USE WITH OPERATING LAPAROSCOPE WITH 6MM Ø INSTRUMENT CHANNEL
LIST NO. 52205-501	6 MM Ø, 33 CM SILASTIC RING APPLICATOR FOR SINGLE AND DUAL RING APPLICATION; USE WITH LAPAROSCOPE AND 6MM Ø TROCAR SLEEVE
LIST NO. 52208-501	8 MM Ø, 27 CM SILASTIC RING APPLICATOR- FOR SINGLE AND DUAL RING APPLICATION; USE WITH LAPAROSCOPE AND 8MM Ø TROCAR SLEEVE SAME DIAMETER AS ACMI
LIST NO. 52208-500	8 MM Ø, 33 CM SILASTIC RING APPLICATOR- FOR SINGLE AND DUAL RING APPLICATION; USE WITH LAPAROSCOPE AND 8MM Ø TROCAR SLEEVE SAME DIAMETER AS ACMI
LIST NO. 52005-501	SILASTIC RING APPLIER AND GUIDE CONE SET [NOTE: ONE (1) SET COMES WITH NEW INSTRUMENT]

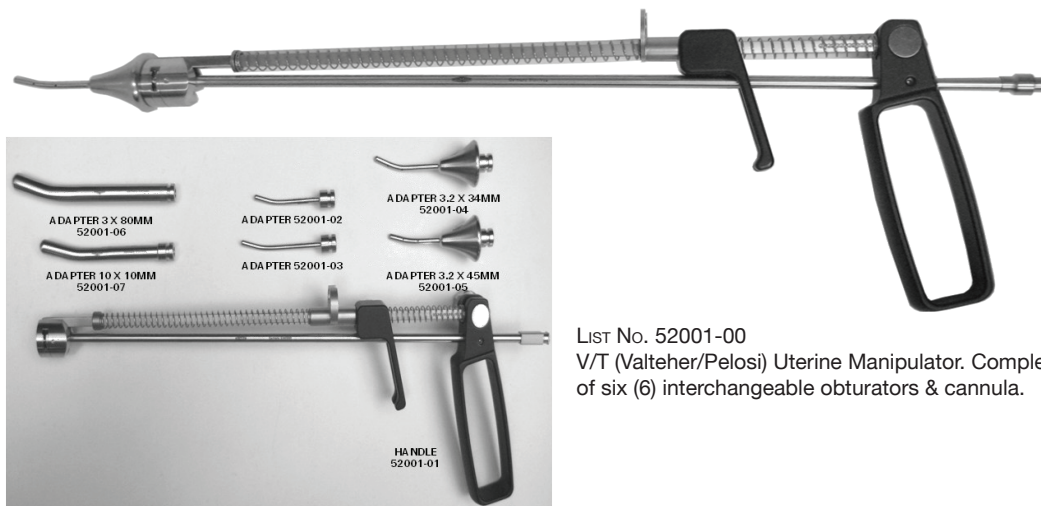

LIST No. 52001-00
V/T (Valteher/Pelosi) Uterine Manipulator. Complete with set of six (6) interchangeable obturators & cannula.



**MILLENNIUM
SURGICAL**

SURGICALINSTRUMENTS.COM

UNMATCHED SERVICE · UNLIMITED INSTRUMENTS

822 Montgomery Ave #205
Narberth, PA 19072
(800) 600-0428


Fallopian Ring Applicator

The Silastic ring is an inert, radio opaque silicone band that is simply applied to a looped section of the fallopian tube with the Silastic Ring Applicator.

The procedure provides high potential for reversibility. The success of surgical anastomosis is most dependent upon the skill of the surgeon.

Illustration of the Tubal Ligation Procedure:

1. Single puncture using 12 mm operating laparoscope.
2. Double puncture requiring a 10 mm laparoscopic and 7 mm Ø trocar sleeve.
3. The Silastic Ring Applicator contains forceps to grasp the fallopian tube 3 to 4 cm from the cornual area.
4. The fallopian tube is gently drawn into the inner cylinder of the instrument, forming a knuckle. The Silastic Ring located on the inner cylinder is released onto the knuckle, occluding the base. This procedure is repeated on the second fallopian tube. Once in place, the Silastic Ring will not slip.

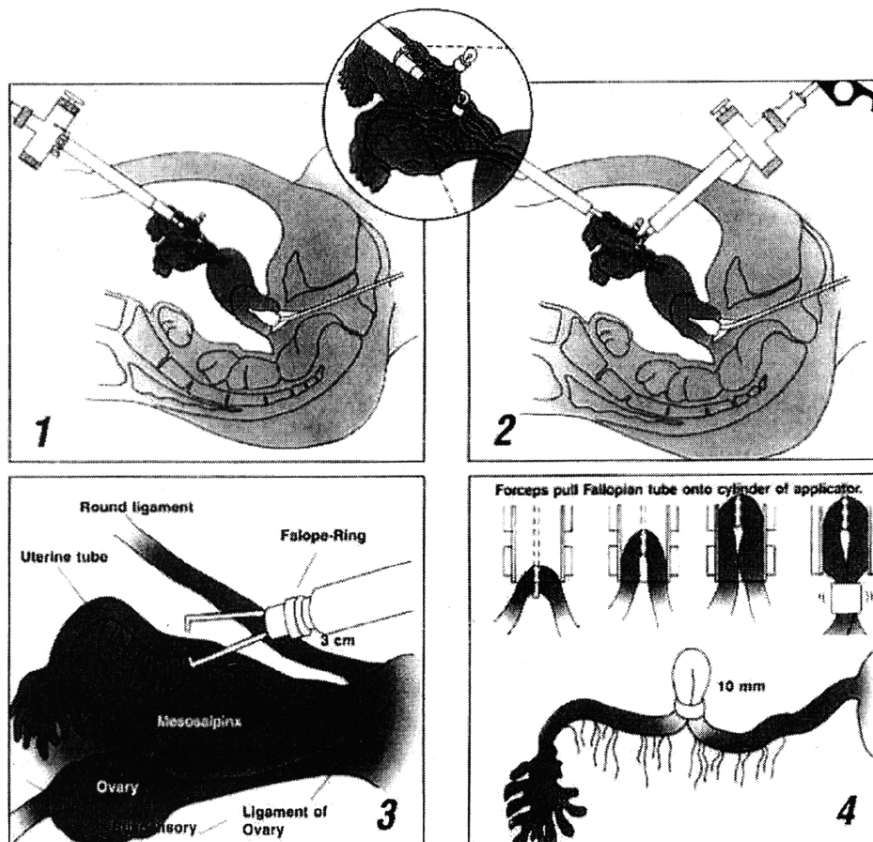

72-52005501

LIST No. 72-52005501

Silastic Ring Applicator and Guide Cone Set

[Note: One (1) set comes with new instrument]